

Rapid Screening of Metals Using Portable High-Resolution X-Ray Fluorescence Spectrometers

Alan D. Hewitt

April 1995



DISTRIBUTION STATEMENT A

Approved for public release; Distribution Unlimited

19950926 105

DTIC QUALITY INCPECTED 5

Abstract

Analysis of copper, zinc, arsenic, lead, chromium, cobalt, nickel, mercury, thallium, selenium, silver, antimony, cadmium, tin, and barium was performed on soils and other particle matrices using two field-portable high-resolution X-ray fluorescence spectrometers (XRF). Quantitative determinations were based on fundamental parameter analysis and a second method that relies on analyte response factors and uses the Compton K_{α} incoherent backscatter peak for matrix normalization. These two methods of instrumental analysis require only a few reference materials and are relatively insensitive to sample matrix composition. This study assessed the capability of these two rapid XRF analysis methods by determining metal concentrations in reference materials, field samples, and laboratory spiked soils. With the exception of nickel, cobalt, and chromium, concentrations within 50% of the expected values were consistently obtained at and below 1000 μ g/g.

For conversion of SI units to non-SI units of measurement consult ASTM Standard E380-93, Standard Practice for Use of the International System of Units, published by the American Society for Testing and Materials, 1916 Race St., Philadelphia, Pa. 19103.

This report is printed on paper that contains a minimum of 50% recycled material.

Special Report 95-14



Rapid Screening of Metals Using Portable High-Resolution X-Ray Fluorescence Spectrometers

Alan D. Hewitt April 1995

Accesio	n For			
NTIS DTIC Unanno Justific	TAB ounced	X		
By Distribu	ution /			
Availability Codes				
Dist		and / or ecial		
A-1				

PREFACE

This report was prepared by Alan D. Hewitt, Research Physical Scientist, Geological Sciences Division, Research and Engineering Directorate, U.S. Army Cold Regions Research and Engineering Laboratory, Hanover, New Hampshire.

Funding for this work was provided by the U.S. Army Environmental Center, Martin H. Stutz, Project Monitor. The author thanks Dr. T.M. Spittler for suggesting the response factor/Compton K_{α} peak normalization method of analysis, Scott Clifford for providing the field samples, and Thomas Ranney and Dr. Thomas F. Jenkins for critical review of the text.

This publication reflects the views of the author and does not suggest or reflect policy, practices, programs, or doctrine of the U.S. Army or of the Government of the United States.

The contents of this report are not to be used for advertising or promotional purposes. Citation of brand names does not constitute an official endorsement or approval of the use of such commercial products.

CONTENTS

		Page
Prefac	e	ii
Introd	uction	1
Instru	mentation	2
	ation	2
	Sis	3
	imental	4
	erence materials	4
	ld samples	4
	ated laboratory soils	4
	s and discussion	4
	usions	11
	ture cited	12
	act	15
Abstra		
ILLUS	STRATIONS	
Figure		
1.	Spectrum of several characteristic K_{α} and K_{β} peaks for metals and	
	Compton and Rayleigh backscatter	3
2.	Lead concentrations established for spiked soils by both RF/Comp. K_{α}	
	normalization and FP analysis	7
3.	Antimony concentrations established for spiked soils by both RF/Comp.	
	K _n normalization and FP analysis	8
4.	Fluorescence yield vs. atomic number for K and L lines	11
5.	Average and standard deviation of concentrations established for	
	soils spiked with 1000 µg metal/g	12
TABL	.ES	
Table		
1.	Primary sources and analyte lines for metals of environmental concern	
	that can be detected by X-ray fluorescence spectrometry	2
2	Daily response factors established for Cu, Zn, As, and Pb based on the	
۷.	SRM 2710 certified reference material	3
3	List of certified reference materials	4
	Characteristics of laboratory-treated soils	4
	Treatment scheme for spiking soil subsamples with metals	Į
	Detection limit estimates and intensity counts	į
	Analysis of commercial reference materials from the National Institute of	
,.	Standards and Technology	(
Q	Analysis of commercial reference materials from Resource Technology	
0.	Corporation	(
Q	Concentrations of Cu, Zn, and Pb determined for field soil samples	
٦.	Constitution of any man and an arrangement of the second o	

Table		Figure
10.	Concentrations of Cr, Cu, Zn, As, and Pb determined for spiked soil	
	matrices by response factor/Compton K_{α} peak normalization	8
11.	Concentrations of Co, Ni, Hg, Tl, and Se determined for spiked soil	
	matrices by response factor/Compton K_{α} peak normalization	9
12.	Concentrations of Ag, Cd, Sn, Sb, and Ba determined for spiked soil	
	matrices by response factor/Compton K_{α} peak normalization	9
13.	Concentrations of Cr, Cu, Zn, As, and Pb determined by fundamental	
	parameter analysis using the Spectrace 9000	10
14.	Concentrations of Co, Ni, Hg, Se, and Fe determined by fundamental	
	parameter analysis using the Spectrace 9000	10
15.	Concentrations of Ag, Cd, Sn, Sb, and Ba determined by fundamental	
	parameter analysis using the Spectrace 9000	11

Rapid Screening of Metals Using Portable High-Resolution X-Ray Fluorescence Spectrometers

ALAN D. HEWITT

INTRODUCTION

Metal pollution over the past decade has come to be of greater concern due to an increasing awareness of the pathways leading to chronic human toxicity (Spittler and Fender 1979). During this time the use of X-ray fluorescence (XRF) analysis has increased for establishing the presence of elevated levels of metals in environmental samples (e.g., Furst et al. 1985, Piorek and Rhodes 1988, Grupp et al. 1989, Watson et al. 1989, Ashe et al. 1991, Carlson and Alexander 1991, Driscoll et al. 1991, Garby 1991, Harding 1991, Puls et al. 1994, etc.). This instrumental method of analysis establishes concentration estimates for several metals over a concentration range extending from percent levels to around $100 \,\mu g/g$ (ppm). The method requires little or no sample pretreatment and can be performed for a small fraction of the time and cost associated with acid extraction/atomic absorption or emission analysis. Moreover, XRF analysis is nondestructive, so subsamples can be archived or analyzed by other procedures. For these reasons, XRF analysis is a very practical and economical method when screening for high concentrations of metals during a remedial investigation and feasibility study (RI/FS) at a suspected hazardous

Today several manufacturers offer energy-dispersive XRF spectrometers that have been designed to be compatible with field operations. The most transportable systems have a battery power supply option and use one or more radioactive sources as the primary incident radiation for elemental excitation. The first generation of these stand-alone systems was equipped with gas proportional detectors that had a spectral resolution on the order of 800–1000 electron volts (1000 eV = 1 keV). More recently, systems have been marketed (HNU Systems Inc., Spectrace Instru-

ments, Inc., Metorex, Inc.) with silicon (drifted with lithium) [Si(Li)] or mercuric iodide (HgI $_2$) detectors capable of achieving resolutions of 170 and 300 eV, respectively. This increased spectral resolution allows for the unambiguous qualitative identification of metals in complex mixtures.

Quantitatively, however, the measurement of discrete spectral energies obtained by XRF analysis is often dependent on other metals present in the sample, due to absorption and enhancement effects. For environmental sample analysis, these matrix-specific effects are best addressed by the use of empirical coefficients (Piorek and Rhodes 1988). This method of instrumental calibration requires the acquisition or development of several well-characterized materials with a matrix composition similar to that of the samples. Moreover, the standards should contain a range of metal concentrations bracketing the desired level of quantitation. Due to these requirements, calibration standards are often both material and site specific, and their preparation and/or verification may take one or more weeks.

A more practical solution for an initial site investigation would be a calibration method that is insensitive to sample matrix so that only few standards would be needed to handle a diverse range of samples (e.g., soil composition "sand/slit/clay," sediment, sludge, dust, paint chips, etc.). Ideally, this capability would allow samples to be analyzed independent of any previous site characterization. With this concept in mind, two quantitation routines have been proposed for field-transportable high-resolution XRF systems. One routine is based on fundamental parameters (FP) analysis, which eliminates the need for several calibration standards by relying on certain physical constants to estimate metal concentration while theoretically correcting for matrix discrepancies (Figura 1987, 1993). Another approach uses response factors (RF)

and corrects for matrix discrepancies by normalizing to the Compton K_{α} (Comp. K_{α}) incoherent radiation backscatter peak.*

This study evaluates these two methods of rapid sample analysis by determining the concentrations of copper (Cu), zinc (Zn), arsenic (As), lead (Pb), chromium (Cr), cobalt (Co), nickel (Ni), mercury (Hg), thallium (Tl), selenium (Se), silver (Ag), antimony (Sb), cadmium (Cd), tin (Sn), and barium (Ba) in a variety of solid-particle matrices. All analyses are performed with transportable highresolution XRF systems that can be configured for stand-alone operations. This combination of rapid sample analysis and instrumental transportability is well suited for screening purposes, thus a data quality objective of an accuracy of $\pm 50\%$ and detection limits of less than $1000~\mu g/g$ was used (Raab et al. 1987).

INSTRUMENTATION

Explanations of the principals of XRF analysis can be found elsewhere (Driscoll et al. 1991, Hewitt 1994a, b). The instruments used in this study were the X-Met 920 (Metorex, Inc.; formerly Outokumpu Electronics) and the Spectrace 9000 (Spectrace Instruments) X-ray spectrophotometers. These two field-portable systems are equipped with surface analysis probes allowing for either in situ or intrusive sample analysis and have software-supported FP analysis capabilities. The Spectrace 9000 is equipped with three primary radioactive sources, Fe-55, Cd-109, and Am-241, and has a HgI₂ solidstate detector. The X-Met 920 has two radioactive sources, Cd-109 and Am-241, and has a Si(Li) detector that requires liquid N_2 for operation. Table 1 lists some of the metals that can be determined by these XRF systems. The table also shows those metals that could be quantitated by FP software programs supplied by the respective manufacturers.

CALIBRATION

The FP software routines used in this study are proprietary to the instrument manufacturer. In general, these FP programs are a series of algorithms based on intensities measured for pure elements, sensitivity coefficients, summations of absorption-enhancement terms in appropriate matrices, and correction factors for overlapping peaks (Figura 1987). The X-Met 920 FP program also required the analysis of a fully characterized standard material, which for this study was SRM 2710, from the National Institute of Science and Technology (NIST).

To perform the RF/Comp. K_{α} normalization analysis, a standard must be available that contains the analytes of interest in a matrix that is physically consistent (dry particles) with the samples. When possible, the RF should be established for well resolved K_{α} , K_{β} , L_{α} , or L_{β} spectral lines (Table 1). Figure 1 is an example of an XRF spectrum obtained with the X-Met 920 showing several characteristic peaks for metals along with the incoherent (Compton) and coherent (Rayleigh) sample matrix backscatter.

For this study a finely ground soil, SRM 2710 certified reference soil from NIST, was used to establish the analyte RFs for Cu, Zn, As, and Pb. Table 2 is an example of some of the daily RFs established for these four metals. For the determination of Cr, the NIST SRM 2711 was spiked with $4000 \, \mu g$

Table 1. Primary sources and analyte lines for metals of environmental concern that can be detected by X-ray fluorescence spectrometry.

		Emission lines (keV)				
Source	Metals	K_{α}	K_{β}	L_{α}	L_{eta}	
FE-55	Crt	5.41	5.95			
Cd-109	Cr†*	5.41	5.95			
Cd-109	Mnt	5.89	6.49			
Cd-109	Fet*	6.40	7.06			
Cd-109	Cot	6.92	7.65			
Cd-109	Nit	7.47	8.30			
Cd-109	Cut*	8.04	8.94			
Cd-109	Znt*	8.63	9.61			
	Ast*	10.5	11.8			
Cd-109	Se†	11.2	12.6			
	Hgt			9.98	11.9	
Cd-109	Τĺ			10.3	12.3	
	Pb†*			10.5	12.6	
Am-241	Ag†	22.1	25.2			
Am-241	Cdt	23.1	26.4			
Am-241	Sn†	25.2	28.8			
Am-241	Sb†	26.2	30.1			
Am-241	Bat	32.0	36.8			

Spectrace 9000 preprogrammed for fundamental parameter analysis of this metal.

^{*} T.M. Spittler, U.S. Environmental Protection Agency, Environmental Services Division, Region 1, Lexington, Massachusetts.

X-Met 920 preprogrammed for fundamental parameter analysis of this metal.

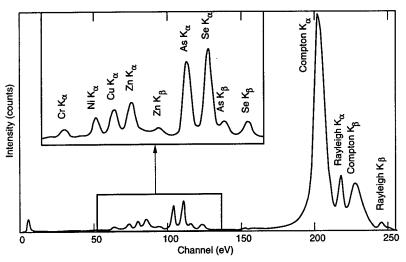


Figure 1. Spectrum of several characteristic K_{α} and K_{β} peaks for metals and Compton and Rayleigh backscatter.

Table 2. Daily response factors established for Cu, Zn, As, and Pb based on the SRM 2710 certified reference material.

Response factor (intensity/concentration)

	Cu	Zn	As	Pb
	12.0	33.9	1.09	39.5
	11.9	34.8	0.90	38.0
	12.3	34.4	1.04	39.6
Average	12.1	34.4	1.01	39.0
Std. dev.	0.21	0.45	0.098	0.90
% RSD	1.7	1.3	9.7	2.3
	SRM	7 2710 conc	entration (µ	g/g)
	2950	6952	626	5532

Cr/g, to establish an RF. All of the other metals analyzed were quantified relative to RFs obtained for the Rocky Mountain Arsenal (RMA) soil matrix spiked with 1000 µg metal/g. Even with a high-resolution XRF instrument, the As K_{α} and Pb L_{α} lines (10.532 and 10.549 keV, respectively) overlap, so the As K_{β} and Pb L_{β} peaks were used. In addition, because high levels (>0.4%) of Fe were present in many of the laboratory-treated soil matrices, the Fe K_{β} peak (7.06 keV) had to be subtracted from the Co K_{α} peak (6.92 keV), due to a spectral overlap. In the case of Co, the K_{β} line was not used because of insufficient response (intensity).

ANALYSIS

Only intrusive sample-analysis procedures were used in this study. Samples consisted of 2- to

5-gram quantities of air-dried materials with particle sizes averaging <600 μm in diameter. All samples were placed into 31-mm-diameter cups and covered with a 0.2-mil polypropylene X-ray film window. XRF analyses were performed for a period of 300 seconds. Field operations using these two preparation steps, i.e., drying and sieving, have shown that up to a hundred samples can be processed and analyzed within a single day (Grupp et al. 1989, Garby 1991).

Fundamental parameter sample analysis was performed using both XRF instruments after the selection of the appropriate software application and assessing

the instrumental tuning. For the Spectrace 9000 this involved checking the resolution and the background response with supplied reference materials and selecting the fine particle application. The FP software application supplied with the X-Met 920 system was based upon the analysis and total characterization of the finely ground NIST SRM 2710.

The RF/Comp. K_{α} normalization analysis method was only performed with the X-Met 920. To perform this method of analysis, the energy spectrum from each analysis was saved and transformed from a 2048- to a 256-channel spectrum for close examination. This approach was used because it was easier to measure the intensity of smoothed peaks and, once transformed, up to six spectra could be overlaid. Measurements included the intensity of elemental spectral lines and the intensity of the incoherent radiation backscatter peak (Comp. K_{α}). The intensity of the incoherent radiation backscatter reflects both the composition of light elements (e.g. carbon, nitrogen, silicon, etc.) in the sample matrix as well as the overall concentration of detectable analytes (Nielson and Sanders 1983). Intensities of both the baseline and peaks of interest were recorded from the computer display after selecting the appropriate spectrum channels. Concentration estimates were then manually determined by multiplying the baseline-corrected analyte signal intensity by the normalization factor, followed by the response factor. The normalization factor is the quotient of the Compton K_{α} peak intensity of the certified reference material, divided by that of the samples. A more complete

description of this analysis has been published elsewhere (Hewitt 1994a, 1994b).

EXPERIMENTAL

Several reference materials of different solid-waste matrices, along with field-contaminated and laboratory-treated soils, were analyzed to assess the performance of both the FP and RF/Comp. K_{α} normalization.

Reference materials

Table 3 lists the certified reference materials purchased from NIST and the Resource Technology Corporation (RTC). Those purchased from NIST have certified concentrations for the total amount of metal present, whereas the RTC materials report certified values based on the USEPA SW846, 3000-series metal acid extraction procedures (U.S. EPA 1986).

Field samples

Six river sediment subsamples from a suspected hazardous-waste site were analyzed. These field subsamples were taken from thoroughly homogenized samples that had been previously characterized by inductively coupled plasma (ICP) analysis following 3000-series metal acid extraction procedures (U.S. EPA 1986).

Treated laboratory soils

Six different soils were spiked with Cr, Cu, Zn, As, Pb, Ni, Se, Hg, Tl, Co, Ag, Cd, Sn, Sb, and Ba by using concentrated 10,000 mg/L aqueous pure element atomic absorption standards (AESAR/

Table 3. List of certified reference materials.

National Institute of Standards and Technology
CDV 4445 Pr
SRM 1645—River sediment
SRM 1646—Estuarine sediment
SRM 1579a—Powdered lead-based paint
SRM 2704—Buffalo River sediment
SRM 2709—San Joaquin soil
SRM 2710—Montana soil
SRM 2711—Montana soil

Resource Technology Corporation

CRM012—Incinerated sludge CRM014—Baghouse dust

CRM013—Paint chips

CRM020—Soil (from EPA Superfund site)

CRM021—Soil (from contaminated waste site)

Table 4. Characteristics of laboratory-treated soils.

Matrix	Sand (%)	Silt and clay (%)	Grain size* (μm)	Wt.† (g)
Ottawa sand	100		400	4
Rocky Mountain Arsenal	NA	NA	NA	4
Lebanon landfill	45	55	300	4
CRREL soil	NA	NA	NA	4
Tampa Bay sediments	95	5	200	4
Ft. Edwards clay	30	70	30	2

- * 95% cut off
- [†] Weight of soil subsample spiked

NA Not analyzed

Alfa, Johnson Matthey). A complete description of how these soil subsamples were treated has been presented elsewhere (Hewitt 1994b). The soil characteristics and the weight of the treated subsamples are shown in Table 4.

Briefly, the soils were air dried and thoroughly mixed prior to placing subsamples into 31-mmdiameter analysis cups. Analyte spikes were made by pipetting between 0.4 and 0.025 mL quantities of the aqueous standards directly onto the individual soil subsamples, increasing the metal concentration by 1000, 500, 250, 125, or 0 μ g/g. Only five analytes were applied to a set of five replicate soil subsamples in order to limit the total volume of solution added to each subsample (0.2 mL/g). Following this protocol, one set of 30 soil subsamples (5×6) was spiked with Cr, Cu, Zn, As, and Pb, a second with Ni, Se, Hg, Tl, and Co, and a third with Ag, Cd, Sn, Sb, and Ba. Analyte additions were performed so that a given soil subsample was not treated with the same concentration more than once (Table 5). A sixth untreated subsample of each soil type was also analyzed with each group and served as the matrix blank.

RESULTS AND DISCUSSION

Table 6 lists the estimates of detection reported by the manufacturers for FP analysis, along with some values established using the RF/Comp. K_{α} normalization and method detection limit (MDL) (Federal Register 1984). This table also includes the analyte intensities measured for the RMA soil treated with 1000 µg metal/g. These intensities were included to provide a means of predicting a detection limit, based on the assumption that there is a fairly constant inverse relationship between these two parameters. This table indicates that regardless of the method of analysis, all of these

Table 5. Treatment scheme for spiking soil subsamples with metals.

Sets	Metal groups						
1	Cr	Cu	Zn	As	Pb		
2	Ni	Se	Hg	Tl	Co		
3	Sb	Ag	Ba	Cd	Sn		
		Treatm	ent concen	trations			
Subsam	ple						
S1	1000*	125	0	500	250		
S2	500	250	1000	0	125		
S3	250	0	125	1000	500		
S4	125	1000	500	250	0		
S5	0	500	250	125	1000		
Matrix							
blank	NF	NF	NF	NF	NF		

^{*} μg/g NF Not fortified

Table 6. Detection limit estimates and intensity counts.

	Detection limits (μg/g)				
Source	Metal	Spectrace 9000*	X-Met 920 [†]	RF/ Compt.**	Peak intensity*†
FE-55	Cr	180	_	_	
Cd-109	Cr	525	325	270	1.18
Cd-109	Mn	410	225		
Cd-109	Fe	225	200		
Cd-109	Co	205	180		1.33
Cd-109	Ni	125	175		3.82
Cd-109	Cu	90	175	54	6.58
Cd-109	Zn	7 0	160	90	6.47
Cd-109	As	50	140	42	3.18
Cd-109	Se	35	140		18.0
Cd-109	Hg	60	_		5.33
Cd-109	ΤĬ			_	7.83
Cd-109	Pb	30		48	8.34
Am-241	Ag	_	7 0	_	25.9
Am-241	Cď	180	100	_	24.6
Am-241	Sn	100	80		27.1
Am-241	Sb	65	80		33.5
Am-241	Ba	20	100	_	31.4

- * Minimum detection limit
- † Minimum determination limit
- ** Method detection limit
- ** Matrix-corrected peak intensity for 1000-ppm spiked RMA soil

metals should be easily quantitated below 1000 μ g/g by XRF analysis.

Tables 7 and 8 show the concentration estimates obtained for the commercial reference materials by these two methods of rapid sample analysis. These two tables show results for Cu, Zn, As, and

Pb when the certified concentration was above one of the estimates of detection listed in Table 6. In only two cases were the values established by the RF/Comp. K_{α} normalization method off by more than ±50% from the certified values. An apparently high concentration was established for Cu in the CRM 021 soil, and a low one for the SRM 2704 river sediment. Fundamental parameter analysis with the Spectrace 9000 XRF failed to establish concentrations within ±50% six times, twice each for Cu, As, and Pb. No values were obtained for Cu in the SRM 2704 river sediment and SRM 2711 soil, while low determinations of As occurred for both the SRM 2710 soil and the CRM 020 soil, and high estimates were obtained for Pb in the CRM 013 paint chips and the CRM 014 baghouse dust. The FP determinations made with the X-Met 920 failed to meet this criterion in only three cases. A low value was established for As in the CRM 020 and high values for Pb in both the CRM 013 paint chips and the CRM 014 baghouse dust.

The high values obtained for the CRM reference materials by these two methods of XRF analysis are not necessarily incorrect, since the certified value is based on an acid extraction that does not necessarily represent the total amount present. However, a low determination for these standards, or one that fails to be within ±50% of the value stated for the NIST reference materials, would be aberrant. The low As concentrations determined by FP analysis had previously been identified as a problem when samples contain much larger (>10 times) quantities of Pb (Harding 1991). The falsenegative Cu determinations were only for samples with certified concentrations very close to the estimates of detection. Overall, these two rapid methods of analysis showed that they were fairly insensitive to this wide variety of particulate matrices by establishing concentrations that would be appropriate for the data quality objectives stated.

Table 9 shows XRF concentration estimates obtained for Cu, Zn, and Pb along with the values obtained by acid-extraction/ICP analysis. The FP analysis was only performed with the X-Met 920 XRF analyzer. In those cases where the values obtained by acid-extraction/ICP analysis were above the appropriate estimate of detection (Table 6), only Cu in sample A and Zn in sample F had XRF estimates that were off by more than $\pm 50\%$. A high Cu value was established by RF/Comp. K_{α} normalization analysis, and a high Zn value was established by FP analysis. Again, since these reported estimates were higher than those obtained after acid extraction and ICP analysis, which is not nec-

essarily a total concentration, the XRF values may not be incorrect.

The concentration estimates for Cu, Zn, As, Pb, Cr, Co, Ni, Hg, Tl, Se, Ag, Cd, Sn, Sb, and Ba in the laboratory-spiked soil subsamples appear in Tables 10 through 15. Tables 10–12 show the values for RF/Comp. K_{α} normalization analysis, and Tables 13–15 show the values for FP analysis using the Spectrace 9000 XRF analyzer. In addition to these metals that were spiked onto the soils, values for iron (Fe) determined by FP analysis are reported in Table 14. Thallium was not determined by FP analysis because this metal had not been included in the software program. All of these concentration estimates were established by using either Cd-109 or Am-241 incident radiation (Table 1).

The values in these tables were corrected for

Table 7. Analysis of commercial reference materials from the National Institute of Standards and Technology.

	M	Metal concentrations (μg/g)			
Standard	Си	Zn	As	Pb	
SRM 1645	109**	1720**		714**	
river sediment	[72]	[1900]		[638]	
	114*	1760*		606*	
SRM 1646		138**			
estuarine sedime	ent	[127]			
SRM 1579a				119,950**	
lead-based paint				[160,000]	
				144,000 [†]	
				124,000*	
SRM 2704	99**	438**		161**	
river sediment	[31]	[302]		[130]	
	ND^{\dagger}	464 [†]		133 [†]	
		527*		200*	
SRM 2709 soil		106**			
		[91]			
		129 [†]			
SRM 2710 soil	2950**	6952**	626**	5532**	
	3350 [†]	7570 [†]	149 [†]	62 4 0 [†]	
SRM 2711 soil	114**	350**	105**	1162**	
	[167]	[343]	[138]	[1100]	
	ND [†]	410 [†]	NDf	1280 [†]	
		414*		1210*	

^{**} Certified concentration

background concentrations present in the soil matrix and for spectral overlap interferences. These corrections were made when both the subsample in the treatment set (Table 5) with no spiked analyte and the subsample of the untreated matrix were determined to have concentrations that exceeded the respective estimates of detection listed in Table 6. However, in cases where there was no estimate of detection for the RF/Comp. K_{α} normalization method of analysis, the lowest value appearing in this table was used. Corrections for Cu, Zn, and Ba were necessary because these metals were present at detectable levels in some of the soil matrices. Corrections were required for Co because of a spectral overlap with the K_{β} peak of Fe and for Ag due to spectral overlap with a peak characteristic of the Am-241 primary radiation source.

Table 8. Analysis of commercial reference materials from the Resource Technology Corporation.

	M	trations (μ	s (μg/g)	
Standard	Cu	Zn	As	Pb
CRM012	3015**	635**		120**
incinerated	[2470]	[342]		[114]
sludge	3590 [†]	744 [†]		126 [†]
sia age	3340*	525*		60*
CRM 013				643**
paint chips				[460]
				2160 [†]
				998*
CRM 014				1914**
baghouse dust				[2080]
O				5360t
				4925*
CRM 020 soil	753**	3021**	397**	5195**
	[687]	[4420]	[429]	[5070]
	521 [†]	3550 ⁺	113 [†]	4580 [†]
	684*	3898*	126*	4950*
CRM 021 soil	5086**	574**		
	[8720]	[549]		
	3060 [†]	408 [†]		

^{**} Certified concentration

^[] Concentration estimate based on response factor/ Compton K_{α} peak normalization with X-Met 920

Concentration estimate based on fundamental parameter analysis with Spectrace 9000

Concentration estimate based on fundamental parameter analysis with X-Met 920

ND Not detected

^[] Concentration estimate based on response factor/ Compton K_{α} peak normalization with X-Met 920

[†] Concentration estimate based on fundamental parameter analysis with Spectrace 9000

^{*} Concentration estimate based on fundamental parameter analysis with X-Met 920

ND Not detected

Table 9. Concentrations (μg/g) of Cu, Zn, and Pb determined for field soil samples.

	Metal concentrations (μg/g)				
Field samples	Cu	Zn	Pb		
Α	114**	1140**	253**		
	[176]	[1260]	[229]		
	286*	1500*	318*		
В	776**	1390**	488**		
	[756]	[1310]	[473]		
	984*	1810*	614*		
С	1860**	261**	2060**		
	[1750]	[182]	[1650]		
	1600*	220*	1600*		
D	3960**	735**	546**		
	[3270]	[408]	[475]		
	3510*	688*	494*		
E	449**	1260**	350**		
	[443]	[1230]	[328]		
	567*	1710*	416*		
F	104**	182**	87**		
	[145]	[200]	[108]		
	331*	391*	245*		

- ** Concentrations based on acid-extraction/ ICP analysis
- [] Concentration estimate based on response factor/Compton K_{α} peak normalization with X-Met 920
 - Concentration estimate based on fundamental parameter analysis with X-Met 920

This study of laboratory-spiked soils was included because commercial reference materials typically contain only a few metals (e.g., Cu, Zn, Pb) at concentrations that can be readily detected by XRF analysis. The soil spiking method used (Hewitt 1994b) appears to have resulted in fairly homogenous and accurate analyte concentrations, as shown by Figures 2 and 3. Plots of the results of Pb and Sb were chosen to show this feature because XRF analysis is particularly sensitive for these two metals (Table 6), and they are representative of metals excited by the Cd-109 and Am-241 primary sources, respectively. It is logical to assume that the other metals spiked onto the different soil matrices were likewise evenly distributed.

The results in Tables 10–15 show that these two methods of rapid sample analysis consistently established concentrations for Cu, Zn, Pb, Se, Ag, Sn, Sb, and Ba that were within $\pm 50\%$ of the expected values from 1000 to 125 ppm. This was also the case for Hg as determined by RF/Comp. K_{α} normalization analysis. The few aberrant values for Hg, as established by FP analysis with the

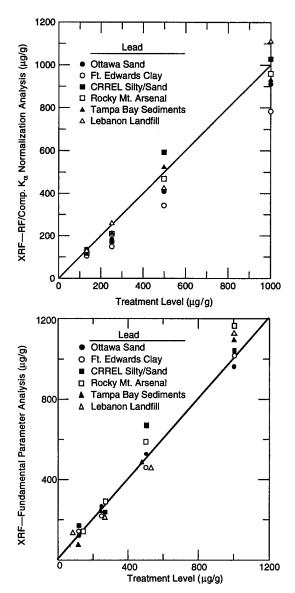


Figure 2. Lead concentrations ($\mu g/g$) established for spiked soils by both RF/Comp. K_{α} normalization and FP analysis.

Spectrace 9000, were most likely caused by a peak overlap with Tl. This problem occurred because Tl had not been one of the metals included in the FP software. With the exception of a few values established for the two lowest treatment levels for As and Cd (and Tl by RF/Comp. K_{α}) these metals were also within the $\pm 50\%$ concentration criterion. Clearly, the concentrations of all of these metals can be adequately estimated at $1000~\mu g/g$ and below by these two rapid methods of XRF analysis. Furthermore, since both methods established concentration trends that were consistent with the treatment levels, they would correctly define the areas of greatest concern.

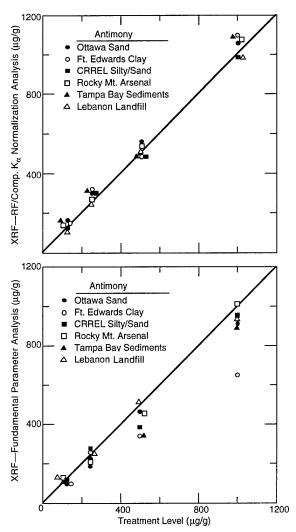


Figure 3. Antimony concentrations ($\mu g/g$) established for spiked soils by both RF/Comp. K_{α} normalization and FP analysis.

The remaining metals in this study, Ni, Co, and Cr, required additional qualifications because their estimates of detection precluded some of the lower laboratory treatment levels. For instance, Ni and Co could only be assessed for spiked concentrations that were >205 ppm $(\mu g/g)$, likewise Cr for concentrations >525 ppm and >270 ppm for FP and the RF/Comp. K_{α} analysis, respectively. Values for all three of these metals were established by both methods of analysis that failed to meet the ±50% concentration criterion. High concentration estimates tended to be established for both Ni and Co by FP analysis in soils with low Fe content (Table 14). In contrast, low concentrations of these two metals were established by RF/Comp. K_{α} normalization analysis when high levels of Fe were present. Both meth-

Table 10. Concentrations (μ g/g) of Cr, Cu, Zn, As, and Pb determined for spiked soil matrices by response factor/Compton K_{α} peak normalization.

	Metal concentrations (μg/g)					
Material & spike	Си	Zn	As	Pb	Cr	
Ottawa sand						
1000 ppm	1010*	1020	1420	911	1770	
500 ppm	534*	458	381	411	856	
250 ppm	264*	224	324	182	312	
125 ppm	132*	127	72	139	84	
0 ppm	69	ND	ND	ND	ND	
Matrix	64	18	ND	10	ND	
Ft. Edwards clay						
1000 ppm	711*	751*	906	785	979	
500 ppm	289*	374*	491	344	638	
250 ppm	144*	252*	229	150	175	
125 ppm	142*	91*	104	108	304	
0 ppm	199	81	ND	19	46	
Matrix	218	80	ND	4	ND	
CRREL soil						
1000 ppm	787*	870	782	1030	1100	
500 ppm	461*	541	635	597	439	
250 ppm	250*	247	320	178	220	
125 ppm	69*	152	26	130	250	
0 ppm	86	48	48	ND	93	
Matrix	62	59	ND	18	257	
Rocky Mountain Ar	senal					
1000 ppm	897*	1050	1450	960	1160	
500 ppm	483*	4 59	608	470	493	
250 ppm	209*	245	304	211	349	
125 ppm	93*	153	105	125	229	
0 ppm	76	ND	3 7	ND	205	
Matrix	62	2	41	23	ND	
Tampa Bay sedimer	nts					
1000 ppm	842*	939	1350*	936	1570	
500 ppm	513*	486	560*	526	653	
250 ppm	250*	218	200*	219	639	
125 ppm	103*	159	73*	124	219	
0 ppm	60	11	220	7	120	
Matrix	71	ND	249	ND	99	
Lebanon landfill so						
1000 ppm	859*	872	988	1192	1500	
500 ppm	373*	469	536	423	495	
250 ppm	195*	224	305	264	407	
125 ppm	94*	140	78	135	326	
0 ppm	99	2	ND	ND	ND	
Matrix	70	8	ND	24	ND	

^{*} Average of 0 ppm and matrix subsamples subtracted

ods of analysis tended to establish high estimates for Cr for soil matrices with low Fe content. Even though this apparent matrix effect was not adequately handled by either of these two methods of analysis to meet the ±50% criterion consistently, the appropriate concentration trends were established.

Table 11. Concentrations ($\mu g/g$) of Co, Ni, Hg, Tl, and Se determined for spiked soil matrices by response factor/ Compton K_{α} peak normalization.

Metal concentrations (µg/g) Material & spike Co Ni Hg TlSe Ottawa sand 1000 ppm 1250 1230 758 909 977 500 ppm 539 609 471 532 502 250 ppm 267 290 149 421 231 199 90 125 ppm 132 225 89 ND ND 52 13 ND 0 ppm ND ND ND 13 Matrix ND Ft. Edwards clay 200* 496 623 813 675 1000 ppm 500 ppm 50* 260 268 482 449 250 ppm ND* 187 145 347 136 ND* 90 120 199 125 ppm 73 6500 45 ND ND 63 0 ppm Matrix 6300 16 ND ND ND CRREL soil 1000 ppm 390* 761 823 993 882 ND* 356 323 553 519 500 ppm 250 ppm 180* 93 141 330 238 125 ppm ND* 15 110 225 128 4050 ND 0 ppm 43 66 8 Matrix 4070 ND ND ND ND Rocky Mountain Arsenal 560* 947 933 1030 977 1000 ppm 500 ppm 270* 481 395 545 483 250 ppm 298 179 60* 449 211 125 ppm 50* 86 161 216 111 0 ppm 1580 12 21 33 ND 1680 ND ND Matrix 12 ND Tampa Bay sediments 1000 ppm 846* 962 933 997 1010 500 ppm 482* 457 395 519 460 178* 315 179 425 250 ppm 221 125 ppm 103* 129 161 199 124 193 22 21 73 ND 0 ppm Matrix 135 ND ND 7 11 Lebanon landfill soil 1000 ppm 435* 838 550 996 968 395* 500 ppm 470 472 460 511 250 ppm 135* 179 207 435 192 85* 125 ppm 88 148 253 123 ND 2260 21 45 ND 0 ppm Matrix 2170 ND ND ND

Several factors that control the sensitivity of elemental analysis by XRF analysis are independent of matrix composition. First of all, elements with low atomic numbers have low fluorescence yields (Fig. 4). The fluorescence yield is the ratio of number of vacancies created within an atom by the incident radiation to the number of vacancies that

Table 12. Concentrations ($\mu g/g$) of Ag, Cd, Sn, Sb, and Ba determined for spiked soil matrices by response factor/Compton K_{α} peak normalization.

	Metal concentrations (μg/g)					
Material & spike	Ag	Cd	Sn	Sb	Ва	
Ottawa sand						
1000 ppm	1170	1000	986	1060	1160	
500 ppm	646	499	398	557	527	
250 ppm	269	176	226	302	262	
125 ppm	109	48	145	160	12	
0 ppm	19	ND	97	ND	NI	
Matrix	25	ND	25	12	NI	
Ft. Edwards clay						
1000 ppm	1032*	936	899	1080	99	
500 ppm	604*	334	447	495	49	
250 ppm	297*	162	248	317	28	
125 ppm	166*	42	146	152	163	
0 ppm	312	ND	92	ND	703	
Matrix	263	ND	26	20	66	
CRREL soil						
1000 ppm	896*	1040	918	996	103	
500 ppm	363*	531	459	493	48	
250 ppm	179*	253	286	295	27	
125 ppm	82*	153	142	142	113	
0 ppm	169	ND	81	ND	39	
Matrix	159	ND	36	14	41	
Rocky Mountain Arse						
1000 ppm	992*	984	745	1070	1060	
500 ppm	410*	416	418	538	530	
250 ppm	297*	142	246	274	220	
125 ppm	140*	42	155	143	130	
0 ppm	246	ND	95	ND	914	
Matrix	230	ND	37	11	96	
Tampa Bay sediments	10101	10/0	000	1000	405	
1000 ppm	1210*	1060	899	1080	1050	
500 ppm	567*	488	447	495	596	
250 ppm	368*	203	348	317	282	
125 ppm	153*	94	146	152	153	
0 ppm	167	ND	92	ND	28	
Matrix	160	ND	26	20	29	
Lebanon landfill soil	1040	1070	003	10/0	004	
1000 ppm	1040	1070	902	1060	990	
500 ppm	535	436	454	496	498	
250 ppm	314	297	248	304	247	
125 ppm	180	182	150	144	118	
0 ppm	49	36	90	ND	360	
Matrix	55	72	31	15	439	

^{*} Average of 0 ppm and matrix subsamples subtracted

actually result in the production of characteristic X-ray photons (Jenkins 1986). In addition, since only a few isotopes are available for field-portable XRF systems, there can be a large separation between incident and excitation energies. The greater this separation, the lower the analyte response, because fewer atoms become excited. A third factor

^{*} Average of 0 ppm and matrix subsamples subtracted

Table 13. Concentrations (μ g/g) of Cr, Cu, Zn, As, and Pb determined by fundamental parameter analysis using the Spectrace 9000.

Metal concentrations (μg/g) PbMaterial & spike Си AsCrZπ Ottawa sand 3010 1000 ppm 1430 1320 1380 976 731 652 389 531 1120 500 ppm 270 287 335 274 580 250 ppm 125 ppm 102 152 51 131 310 ND ND ND ND ND 0 ppm ND Matrix ND ND 10 ND Ft. Edwards clay 1140* 1020 850 1340 1020 1000 ppm 524* 660 467 300 532 500 ppm 334 192* 343 225 ND 250 ppm 147* 125 ppm 227 13 138 ND 0 ppm 131 99 ND ND ND Matrix 67 99 66 ND ND CRREL soil 1280 1040 940 1000 ppm 1130 1270* 57 489 550* 656 657 500 ppm 253 221* 311 232 53 250 ppm 172 77 120* 60 ND 125 ppm 0 ppm ND 74 ND ND ND Matrix ND 95 56 ND ND Rocky Mountain Arsenal 1280 1240 1350 1170 1480 1000 ppm 500 ppm 678 675 636 589 440 333 375 281 250 ppm 236 280 74 125 19 141 120 125 ppm 0 ppm ND ND 22 ND ND 28 ND ND 36 ND Matrix Tampa Bay sediments 1210 1250 1230 1100 1490 1000 ppm 651 477 700 500 ppm 662 545 256 250 ppm 312 289 278 120 83 ND 125 ppm 70 118 63 ND ND 12 ND 66 0 ppm ND Matrix ND ND ND Lebanon landfill soil 1000 ppm 1130 1130 1110 1130 1350 584 567 450 613 460 500 ppm 250 ppm 316 327 332 218 150 125 ppm 71 133 59 130 ND 0 ppm . . ND 13 11 ND ND

Matrix

ND

22

ND

ND

is the depth of penetration achieved by the incident radiation, i.e., how much of the sample is actually analyzed. The penetration depth of incident radiation is inversely proportional to its wavelength and directly proportional to its energy. For the sources used in this study, the depth of penetration (99% attenuation) for a quartz matrix

Table 14. Concentrations (μ g/g) of Co, Ni, Hg, Se, and Fe determined by fundamental parameter analysis using the Spectrace 9000.

		Metal co			<u> </u>
Material & spike	Со	Ni	Hg	Se	Fe
Ottawa sand					
1000 ppm	1830	1538	1010	1240	NI
500 ppm	994	846	781	637	7
250 ppm	442	427	228	307	6
125 ppm	109	183	380	143	12
0 ppm	ND	ND	111	ND	11
Matrix	ND	ND	ND	ND	142
Ft. Edwards clay					
1000 ppm	727*	999	791	1100	55,40
500 ppm	277*	495	458	503	54,60
250 ppm	67*	207	200	265	56,70
125 ppm	ND*	12	360	120	57,00
0 ppm	560	ND	80	ND	55,60
Matrix	406	ND	ND	ND	56,90
CRREL soil					
1000 ppm	965*	1210	1321	1350	40,00
500 ppm	295*	624	665	611	40,90
250 ppm	235*	216	259	332	41,10
125 ppm	315*	148	155	146	41,10
0 ppm	190	11	95	ND	41,30
Matrix	280	ND	16	ND	40,90
Rocky Mountain Ar	senal				
1000 ppm	1490	1480	1170	1390	18,60
500 ppm	900	786	736	689	19,60
250 ppm	479	324	301	304	22,50
125 ppm	224	141	579	140	21,70
0 ppm	ND	ND	121	ND	21,10
Matrix	ND	ND	ND	ND	20,10
Tampa Bay sedimen	ıts				
1000 ppm	1 7 10	1720	1350	1310	441
500 ppm	738	687	728	758	424
250 ppm	423	308	356	350	393
125 ppm	198	152	529	185	347
0 ppm	122	32	151	ND	350
Matrix	24	ND	ND	ND	419
Lebanon landfill soi	1				
1000 ppm	1330	1280	1040	1260	24,90
500 ppm	660	578	744	705	25,90
250 ppm	208	294	227	346	25,50
125 ppm	160	119	336	134	27,20
0 ppm	ND	ND	151	ND	25,60
Matrix	109	14	ND	ND	25,70

^{*} Average of 0 ppm and matrix subsamples subtracted

ranges from 40 mm for Am-241 to 4.8 mm for Cd-109. The combination of these factors controls the intensity that will be measured for a given element. For example, Table 6 shows the range of analyte intensities obtained for the same concentration in the RMA soil.

To illustrate the combination of these two ef-

^{*} Average of 0 ppm and matrix subsamples subtracted

Table 15. Concentrations (μ g/g) of Ag, Cd, Sn, Sb, and Ba determined by fundamental parameter analysis using the Spectrace 9000.

Material & spike Ag Cd Sn Sb Ba Ottawa sand 1000 ppm 937 1110 1000 915 772 500 ppm 504 581 429 474 353 250 ppm 274 394 269 196 174 125 ppm 193 229 174 109 75 0 ppm 77 111 43 ND ND Matrix 25 74 41 ND ND Ft. Edwards clay 1000 ppm 820 976 857 654 601* 500 ppm 377 542 452 345 300* 250 ppm 197 361 226 261 215* 125 ppm 197 361 226 261 215* 125 ppm 197 361 226 261 215* 125 ppm ND 100 150 955 715* 500 ppm <t< th=""><th></th><th colspan="5">Metal concentrations (μg/g,</th></t<>		Metal concentrations (μg/g,				
1000 ppm	Material & spike	Ag	Cd	Sn	Sb	Ва
1000 ppm	Ottowa cand					
500 ppm 504 581 429 474 353 250 ppm 274 394 269 196 174 125 ppm 193 229 174 109 75 0 ppm 77 111 43 ND ND Matrix 25 74 41 ND ND Ft. Edwards clay 1000 ppm 820 976 857 654 601* 500 ppm 377 542 452 345 300* 250 ppm 197 361 226 261 215* 125 ppm 86 172 154 108 148* 0 ppm ND 127 74 ND 379 Matrix ND 100 81 17 436 CRREL soil 1000 ppm 998 1060 1050 955 715* 500 ppm 508 584 527 490 381* 225 ppm 201* 25		937	1110	1000	915	772
250 ppm	* *					
125 ppm 193 229 174 109 75 0 ppm 77 111 43 ND ND ND Matrix 25 74 41 ND ND ND Matrix 25 74 41 ND ND ND Matrix 25 74 41 ND ND ND ND ND Matrix 25 74 41 ND						
O ppm	* *					
Matrix 25 74 41 ND ND Ft. Edwards clay 1000 ppm 820 976 857 654 601* 500 ppm 377 542 452 345 300* 250 ppm 197 361 226 261 215* 125 ppm 86 172 154 108 148* 0 ppm ND 127 74 ND 379 Matrix ND 100 81 17 436 CRREL soil 1000 ppm 998 1060 1050 955 715* 500 ppm 508 584 527 490 381* 250 ppm 268 302 301 276 201* 125 ppm 191 225 153 122 92* 0p* 9p* 381* 227 290 381* 250 pp* 9p* 486 562 472* 29* 0p* 381* 276 201* 29*						
1000 ppm 820 976 857 654 601* 500 ppm 377 542 452 345 300* 250 ppm 197 361 226 261 215* 125 ppm 86 172 154 108 148* 0 ppm ND 127 74 ND 379 Matrix ND 100 81 17 436 CRREL soil 1000 ppm 998 1060 1050 955 715* 500 ppm 508 584 527 490 381* 250 ppm 268 302 301 276 201* 125 ppm 191 225 153 122 92* 0p* 0p* 388 Matrix 40 99 28 ND 329 Rocky Mountain Arsenal 1000 ppm 1060 1200 1040 1010 850* 500 ppm 535 613 486 562						
1000 ppm 820 976 857 654 601* 500 ppm 377 542 452 345 300* 250 ppm 197 361 226 261 215* 125 ppm 86 172 154 108 148* 0 ppm ND 127 74 ND 379 Matrix ND 100 81 17 436 CRREL soil 1000 ppm 998 1060 1050 955 715* 500 ppm 508 584 527 490 381* 250 ppm 268 302 301 276 201* 125 ppm 191 225 153 122 92* 0p* 0p* 388 Matrix 40 99 28 ND 329 Rocky Mountain Arsenal 1000 ppm 1060 1200 1040 1010 850* 500 ppm 535 613 486 562	Ft. Edwards clay					
500 ppm 377 542 452 345 300* 250 ppm 197 361 226 261 215* 125 ppm 86 172 154 108 148* 0 ppm ND 127 74 ND 379 Matrix ND 100 81 17 436 CRREL soil 301 100 81 17 436 CRREL soil 302 301 276 201* 500 ppm 508 584 527 490 381* 250 ppm 268 302 301 276 201* 125 ppm 191 225 153 122 92* 0 ppm 67 51 37 16 288 Matrix 40 99 28 ND 329 Rocky Mountain Arsenal 3100 1040 1010 850* 500 ppm 535 613 486 562 472* 250 ppm 296 406 336 221 250* <tr< td=""><td></td><td>820</td><td>976</td><td>857</td><td>654</td><td>601*</td></tr<>		820	976	857	654	601*
250 ppm		377	542	452	345	300*
125 ppm	• •	197	361	226	261	215*
0 ppm Matrix ND 100 127 74 ND 379 Matrix ND 100 81 17 436 CRREL soil 1000 ppm 998 1060 1050 955 715* 500 ppm 508 584 527 490 381* 250 ppm 268 302 301 276 201* 125 ppm 191 225 153 122 92* 0 ppm 67 51 37 16 288 Matrix 40 99 28 ND 329 Rocky Mountain Arsenal 1000 ppm 1060 1200 1040 1010 850* 500 ppm 535 613 486 562 472* 250 ppm 296 406 336 221 250* 125 ppm 190 253 194 133 105* 0 ppm 77 143 33 22 716 Matrix 37 80 58 14 734 Tampa Bay sediments 1000 ppm 989 1141 1030 897 706 500 ppm 486 557 485 445 365 250 ppm 268 299 205 235 215 125 ppm 181 189 114 115 92 0 ppm 56 42 51 ND ND Matrix 56 84 55 12 ND ND Matrix 56 84 55 12 ND ND 1000 ppm 56 42 51 ND ND ND Matrix 56 84 55 12 ND ND	• •	86	172	154	108	148*
Matrix ND 100 81 17 436 CRREL soil 1000 ppm 998 1060 1050 955 715* 500 ppm 508 584 527 490 381* 250 ppm 268 302 301 276 201* 125 ppm 191 225 153 122 92* 0 ppm 67 51 37 16 288 Matrix 40 99 28 ND 329 Rocky Mountain Arsenal 30 1040 1010 850* 500 ppm 535 613 486 562 472* 250 ppm 296 406 336 221 250* 125 ppm 190 253 194 133 105* 0 ppm 77 143 33 22 716 Matrix 37 80 58 14 734 Tampa Bay sediments 1000 ppm 989		ND	127	74	ND	379
1000 ppm 998 1060 1050 955 715* 500 ppm 508 584 527 490 381* 250 ppm 268 302 301 276 201* 125 ppm 191 225 153 122 92* 0 ppm 67 51 37 16 288 Matrix 40 99 28 ND 329 Rocky Mountain Arsenal 1000 ppm 1060 1200 1040 1010 850* 500 ppm 535 613 486 562 472* 250 ppm 296 406 336 221 250* 125 ppm 190 253 194 133 105* 0 ppm 77 143 33 22 716 Matrix 37 80 58 14 734 Tampa Bay sediments 1000 ppm 989 1141 1030 897	* *	ND	100	81	17	436
500 ppm 508 584 527 490 381* 250 ppm 268 302 301 276 201* 125 ppm 191 225 153 122 92* 0 ppm 67 51 37 16 288 Matrix 40 99 28 ND 329 Rocky Mountain Arsenal 30 28 ND 329 Rocky Mountain Arsenal 30 1040 1010 850* 500 ppm 1060 1200 1040 1010 850* 500 ppm 535 613 486 562 472* 250 ppm 296 406 336 221 250* 125 ppm 190 253 194 133 105* 0 ppm 77 143 33 22 716 Matrix 37 80 58 14 734 Tampa Bay sediments 1000 ppm 989 1141 1030	CRREL soil					
250 ppm 268 302 301 276 201* 125 ppm 191 225 153 122 92* 0 ppm 67 51 37 16 288 Matrix 40 99 28 ND 329 Rocky Mountain Arsenal 1000 ppm 1060 1200 1040 1010 850* 500 ppm 535 613 486 562 472* 250 ppm 296 406 336 221 250* 125 ppm 190 253 194 133 105* 0 ppm 77 143 33 22 716 Matrix 37 80 58 14 734 Tampa Bay sediments 1000 ppm 989 1141 1030 897 706 500 ppm 486 557 485 445 365 250 ppm 268 299 205 235	1000 ppm	998	1060	1050	955	715*
125 ppm 191 225 153 122 92* 0 ppm 67 51 37 16 288 Matrix 40 99 28 ND 329 Rocky Mountain Arsenal 1000 ppm 1060 1200 1040 1010 850* 500 ppm 535 613 486 562 472* 250 ppm 296 406 336 221 250* 125 ppm 190 253 194 133 105* 0 ppm 77 143 33 22 716 Matrix 37 80 58 14 734 Tampa Bay sediments 1000 ppm 989 1141 1030 897 706 500 ppm 486 557 485 445 365 250 ppm 268 299 205 235 215 125 ppm 181 189 114 115 92 0 ppm 56 42 51 ND ND Matrix 56 84 55 12 ND Lebanon landfill soil 1000 ppm 914 1190 1020 930 832* 500 ppm 525 608 538 515 343* 250 ppm 301 337 287 259 232* 125 ppm 301 337 287 259 232* 125 ppm 205 223 165 133 111* 0 ppm 39 66 41 ND 278	500 ppm	508	584	527	490	381*
0 ppm 67 51 37 16 288 Matrix 40 99 28 ND 329 Rocky Mountain Arsenal 1000 ppm 1060 1200 1040 1010 850* 500 ppm 535 613 486 562 472* 250 ppm 296 406 336 221 250* 125 ppm 190 253 194 133 105* 0 ppm 77 143 33 22 716 Matrix 37 80 58 14 734 Tampa Bay sediments 1000 ppm 989 1141 1030 897 706 500 ppm 486 557 485 445 365 250 ppm 268 299 205 235 215 125 ppm 181 189 114 115 92 0 ppm 56 42 51 ND ND<		268	302	301	276	201*
Matrix 40 99 28 ND 329 Rocky Mountain Arsenal 1000 ppm 1060 1200 1040 1010 850* 500 ppm 535 613 486 562 472* 250 ppm 296 406 336 221 250* 125 ppm 190 253 194 133 105* 0 ppm 77 143 33 22 716 Matrix 37 80 58 14 734 Tampa Bay sediments 1000 ppm 989 1141 1030 897 706 500 ppm 486 557 485 445 365 250 ppm 268 299 205 235 215 125 ppm 181 189 114 115 92 0 ppm 56 42 51 ND ND Matrix 56 84 55 12 ND Lebanon landfill soil	125 ppm	191	225	153	122	92*
Rocky Mountain Arsenal 1000 ppm 1060 1200 1040 1010 850* 500 ppm 535 613 486 562 472* 250 ppm 296 406 336 221 250* 125 ppm 190 253 194 133 105* 0 ppm 77 143 33 22 716 Matrix 37 80 58 14 734 Tampa Bay sediments 1000 ppm 989 1141 1030 897 706 500 ppm 486 557 485 445 365 250 ppm 268 299 205 235 215 125 ppm 181 189 114 115 92 0 ppm 56 42 51 ND ND Matrix 56 84 55 12 ND Lebanon landfill soil 11000 ppm 914 1190 1020 930 832* 500 ppm 525 608 538	0 ppm	67	51	37	16	288
1000 ppm 1060 1200 1040 1010 850* 500 ppm 535 613 486 562 472* 250 ppm 296 406 336 221 250* 125 ppm 190 253 194 133 105* 0 ppm 77 143 33 22 716 Matrix 37 80 58 14 734 Tampa Bay sediments 1000 ppm 989 1141 1030 897 706 500 ppm 486 557 485 445 365 250 ppm 268 299 205 235 215 125 ppm 181 189 114 115 92 0 ppm 56 42 51 ND ND Matrix 56 84 55 12 ND Lebanon landfill soil 11000 ppm 914 1190 1020 930 832* <tr< td=""><td>Matrix</td><td>40</td><td>99</td><td>28</td><td>ND</td><td>329</td></tr<>	Matrix	40	99	28	ND	329
500 ppm 535 613 486 562 472* 250 ppm 296 406 336 221 250* 125 ppm 190 253 194 133 105* 0 ppm 77 143 33 22 716 Matrix 37 80 58 14 734 Tampa Bay sediments 1000 ppm 989 1141 1030 897 706 500 ppm 486 557 485 445 365 250 ppm 268 299 205 235 215 125 ppm 181 189 114 115 92 0 ppm 56 42 51 ND ND Matrix 56 84 55 12 ND Lebanon landfill soil 11000 ppm 914 1190 1020 930 832* 500 ppm 525 608 538 515 343*	Rocky Mountain Arsena					
250 ppm 296 406 336 221 250* 125 ppm 190 253 194 133 105* 0 ppm 77 143 33 22 716 Matrix 37 80 58 14 734 Tampa Bay sediments 1000 ppm 989 1141 1030 897 706 500 ppm 486 557 485 445 365 250 ppm 268 299 205 235 215 125 ppm 181 189 114 115 92 0 ppm 56 42 51 ND ND Matrix 56 84 55 12 ND Lebanon landfill soil 1000 ppm 914 1190 1020 930 832* 500 ppm 525 608 538 515 343* 250 ppm 301 337 287 259 232* 125 ppm 205 223 165 133 111* 0 ppm 39 66 41 ND 278						
125 ppm 190 253 194 133 105* 0 ppm 77 143 33 22 716 Matrix 37 80 58 14 734 Tampa Bay sediments 1000 ppm 989 1141 1030 897 706 500 ppm 486 557 485 445 365 250 ppm 268 299 205 235 215 125 ppm 181 189 114 115 92 0 ppm 56 42 51 ND ND Matrix 56 84 55 12 ND Lebanon landfill soil 11000 ppm 914 1190 1020 930 832* 500 ppm 525 608 538 515 343* 250 ppm 301 337 287 259 232* 125 ppm 205 223 165 133 111* 0 ppm 39 66 41 ND 278						
0 ppm 77 143 33 22 716 Matrix 37 80 58 14 734 Tampa Bay sediments 1000 ppm 989 1141 1030 897 706 500 ppm 486 557 485 445 365 250 ppm 268 299 205 235 215 125 ppm 181 189 114 115 92 0 ppm 56 42 51 ND ND Matrix 56 84 55 12 ND Lebanon landfill soil 1000 ppm 914 1190 1020 930 832* 500 ppm 525 608 538 515 343* 250 ppm 301 337 287 259 232* 125 ppm 205 223 165 133 111* 0 ppm 39 66 41 ND 278						
Matrix 37 80 58 14 734 Tampa Bay sediments 1000 ppm 989 1141 1030 897 706 500 ppm 486 557 485 445 365 250 ppm 268 299 205 235 215 125 ppm 181 189 114 115 92 0 ppm 56 42 51 ND ND Matrix 56 84 55 12 ND Lebanon landfill soil 1190 1020 930 832* 500 ppm 914 1190 1020 930 832* 500 ppm 525 608 538 515 343* 250 ppm 301 337 287 259 232* 125 ppm 205 223 165 133 111* 0 ppm 39 66 41 ND 278	125 ppm					
Tampa Bay sediments 1000 ppm 989 1141 1030 897 706 500 ppm 486 557 485 445 365 250 ppm 268 299 205 235 215 125 ppm 181 189 114 115 92 0 ppm 56 42 51 ND ND Matrix 56 84 55 12 ND Lebanon landfill soil 1000 ppm 914 1190 1020 930 832* 500 ppm 525 608 538 515 343* 250 ppm 301 337 287 259 232* 125 ppm 205 223 165 133 111* 0 ppm 39 66 41 ND 278	* *					
1000 ppm 989 1141 1030 897 706 500 ppm 486 557 485 445 365 250 ppm 268 299 205 235 215 125 ppm 181 189 114 115 92 0 ppm 56 42 51 ND ND Matrix 56 84 55 12 ND Lebanon landfill soil 1000 ppm 914 1190 1020 930 832* 500 ppm 525 608 538 515 343* 250 ppm 301 337 287 259 232* 125 ppm 205 223 165 133 111* 0 ppm 39 66 41 ND 278	Matrix	37	80	58	14	734
500 ppm 486 557 485 445 365 250 ppm 268 299 205 235 215 125 ppm 181 189 114 115 92 0 ppm 56 42 51 ND ND Matrix 56 84 55 12 ND Lebanon landfill soil 1000 ppm 914 1190 1020 930 832* 500 ppm 525 608 538 515 343* 250 ppm 301 337 287 259 232* 125 ppm 205 223 165 133 111* 0 ppm 39 66 41 ND 278						
250 ppm 268 299 205 235 215 125 ppm 181 189 114 115 92 0 ppm 56 42 51 ND ND Matrix 56 84 55 12 ND Lebanon landfill soil 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 300 <td< td=""><td>• •</td><td></td><td></td><td></td><td></td><td></td></td<>	• •					
125 ppm 181 189 114 115 92 0 ppm 56 42 51 ND ND Matrix 56 84 55 12 ND Lebanon landfill soil 84 1000 ppm 930 832* 500 ppm 525 608 538 515 343* 250 ppm 301 337 287 259 232* 125 ppm 205 223 165 133 111* 0 ppm 39 66 41 ND 278						
0 ppm Matrix 56 56 42 84 51 55 ND 12 ND ND Lebanon landfill soil 84 55 12 ND 1000 ppm 914 1190 1020 930 832* 500 ppm 525 608 538 515 343* 250 ppm 301 337 287 259 232* 125 ppm 205 223 165 133 111* 0 ppm 39 66 41 ND 278						
Matrix 56 84 55 12 ND Lebanon landfill soil 1000 ppm 914 1190 1020 930 832* 500 ppm 525 608 538 515 343* 250 ppm 301 337 287 259 232* 125 ppm 205 223 165 133 111* 0 ppm 39 66 41 ND 278	* *					
1000 ppm 914 1190 1020 930 832* 500 ppm 525 608 538 515 343* 250 ppm 301 337 287 259 232* 125 ppm 205 223 165 133 111* 0 ppm 39 66 41 ND 278		• •				
1000 ppm 914 1190 1020 930 832* 500 ppm 525 608 538 515 343* 250 ppm 301 337 287 259 232* 125 ppm 205 223 165 133 111* 0 ppm 39 66 41 ND 278	Lebanon landfill soil					
500 ppm 525 608 538 515 343* 250 ppm 301 337 287 259 232* 125 ppm 205 223 165 133 111* 0 ppm 39 66 41 ND 278		914	1190	1020	930	832*
250 ppm 301 337 287 259 232* 125 ppm 205 223 165 133 111* 0 ppm 39 66 41 ND 278						
125 ppm 205 223 165 133 111* 0 ppm 39 66 41 ND 278	• •					
0 ppm 39 66 41 ND 278						
· F F · · ·	• •					
	Matrix	ND	18	71	18	300

^{*} Average of 0 ppm and matrix subsamples subtracted

fects, matrix composition and characteristic analyte response, Figure 5 shows the average concentration and standard deviations for the 1000 µg metal/g spiked soil matrices. Overall, both the precision and accuracy of analysis improves with atomic number. However, for Cr, Co, and Ni, the degree of bias and range of uncertainty shows that values

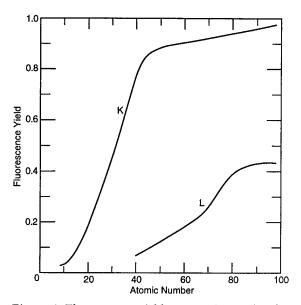


Figure 4. Fluorescence yield vs. atomic number for K and L lines.

even at this level can be greater than ±50% from the true concentration when analyzing a variety of soils. Most likely, the bias in the determined concentrations is the result of a matrix effect, while the extent of uncertainty is a function of both the matrix effect and the characteristic analyte response.

In summary, these two methods of rapid sample analysis with field-portable XRF systems often achieved the data requirement for screening of ±50% accuracy at and below a concentration of 1000 μg/g. This goal was accomplished for a variety of soil matrices, and the results in Tables 7–9 show that this approach is also promising for several other particulate matrices (e.g., sediment, dust, paint chips, and sludge). The determination of Ni, Co, and Cr, however, was not found to accomplish this goal consistently, even though reported estimates of detection for XRF analysis are well below 1000 μg/g. It appears that for analytes with poor XRF sensitivity, matrix effects are more problematic, and matrixspecific standards would be necessary to achieve a ±50% accuracy.

CONCLUSIONS

With the exception of Ni, Co, and Cr, the determination of Cu, Zn, As, Pb, Hg, Tl, Se, Ag, Cd, Sn, Sb, and Ba at and below 1000 μ g/g in a variety of solid particulate materials was often within $\pm 50\%$ of the expected values when using either funda-

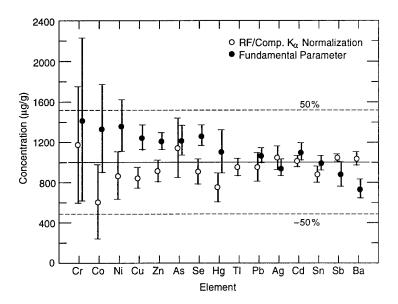


Figure 5. Average and standard deviation of concentrations established for soils spiked with 1000 μg metal/g, as determined by RF/Comp. K_{α} normalization and FP analysis.

mental parameter or response factor/Compton K_{α} peak-normalization methods of analysis. These alternative approaches to XRF analysis are very useful for screening a variety of matrices during RI/FS activities when it is impractical to produce matrix-matched standards.

LITERATURE CITED

Ashe, J.B., P.F. Berry, G.R. Voots, M. Bernick and G. Prince (1991) A high resolution portable XRF HgI₂ spectrometer for field screening of hazardous waste. In *Proceedings of the 2nd International Symposium on Field Screening for Hazardous Waste and Toxic Chemicals*, February 12–14, Las Vegas, Nevada, p. 507–514.

Carlson, C.D. and J.R. Alexander (1991) Data quality assurance/quality control for field X-ray fluorescence spectrometry. In *Proceedings of the 2nd International Symposium on Field Screening for Hazardous Waste and Toxic Chemicals, February 12–14, Las Vegas, Nevada*, p. 525–532.

Driscoll, J.N., J.K. Marshall, C. Wood and T. Spittler (1991) A multifunctional portable X-ray fluorescence instrument. *American Laboratory*, July, p. 25–36.

Federal Register (1984) Definition and procedure for the determination of the method detection limit. Code of Federal Regulations, Part 136, Appendix B, October 26.

Figura, P.M. (1993) Standardless X-ray fluorescence analysis of liquids. *American Laboratory*, July, p. 40–43.

Figura, P.M. (1987) Standardless quantitative X-

ray fluorescence analysis using stored calibration constants. *American Laboratory*, February, **19**(2): 156–164.

Furst, G.A., V. Tillinghast and T.M. Spittler (1985) Screening for metals at hazardous waste sites: A rapid cost-effective technique using X-ray fluorescence. In *Proceedings of the National Conference on Management of Uncontrolled Hazardous Waste Sites, Washington, D.C.*,

Garby, J.C. (1991) Comparison of mobile laboratory XRF and CLP split sample lead results from a Superfund site remediation in New Jersey. In *Proceedings of the 2nd International Symposium on Field Screening for Hazardous Waste and Toxic Chemicals, February 12–14, Las Vegas, Nevada*, p. 671–672.

Grupp, D.J., D.A. Everitt, R.J. Bath, and R. Spear (1989) Use of a transportable XRF spectrometer for on-site analysis of Hg in soils. *American Laboratory*, November, p. 32–40.

Harding, A.R. (1991) Low concentration soil contaminant characterization using EDXRF analysis. In *Proceedings of the 2nd International Symposium on Field Screening for Hazardous Waste and Toxic Chemicals, February* 12–14, *Las Vegas, Nevada*, p. 517–523. **Hewitt, A.D.** (1994a) Screening for metals by X-ray fluorescence spectrometry/response factor/Compton $K_α$ peak normalization analysis. *American Environmental Laboratory*, June, p. 24–26.

Hewitt, A.D. (1994b) Screening for metals by X-ray fluorescence spectrometry using a single calibration standard. USA Cold Regions Research and Engineering Laboratory, Special Report 94-20.

Jenkins, R. (1984) X-ray fluorescence analysis. *Analytical Chemistry*, **56**(9): 1099A–1103A

Nielson, K.K., and R.W. Sanders (1983) Multi-

element analysis of unweighed biological and geological samples using backscatter and fundamental parameters. *Advances in X-Ray Analysis* 1983, **26**:385–390.

Piorek, S. and J.R. Rhodes (1988) A new calibration technique for X-ray analyzers used in hazardous waste screening. In *Proceedings of the 5th National Conference on Hazardous Waste and Hazardous Materials, April 19–21, Las Vegas, Nevada,* p. 428–433. Puls, R., D.A. Clark, C. Carlson and J. Vardy (1994) Characterization of chromium-contaminated soils using field-portable X-ray fluorescence. *Ground Water Monitoring and Remediation*, 14: 111–115.

Raab, G.A., D. Cardenas, S.J. Simon and L.A. Eccles (1987) Evaluation of a prototype field-portable X-ray fluorescence system for hazardous

waste screening. Environmental Monitoring Systems Laboratory, Las Vegas, Nevada, EPA/600/4-87/021, National Technical Information Service, PB87-227633.

Spittler, T.M. and **W.A.** Fender (1979) A study of soil contamination and plant lead uptake in Boston urban gardens. *Communications in Soil Science and Plant Analysis*, **10**(9): 1195–1210.

U.S. EPA (1986) *Test Methods for Evaluating Solid Waste*, Vol. 1B. Washington, D.C.: Environmental Protection Agency.

Watson, W., J.P. Walsh and B. Glynn (1989) Onsite X-ray fluorescence spectrometry mapping of metal contaminants in soils at Superfund sites. *American Laboratory*, July.

REPORT DOCUMENTATION PAGE

Form Approved OMB No. 0704-0188

Public reporting burden for this collection of information is estimated to average 1 hour per response, including the time for reviewing instructions, searching existing data sources, gathering and maintaining the data needed, and completing and reviewing the collection of information. Send comments regarding this burden estimate or any other aspect of this collection of information, including suggestion for reducing this burden, to Washington Headquarters Services, Directorate for Information Operations and Reports, 1215 Jefferson Davis Highway, Suite 1204, Arlington, VA 22202-4302, and to the Office of Management and Budget, Paperwork Reduction Project (0704-0188), Washington, DC 20503.

1. AGENCY USE ONLY (Leave blank)	2. REPO	RT DATE April 1995	3. REPORT TY	PE AND DA	TES COVERED		
4. TITLE AND SUBTITLE		71pm 1995		5. FUNDI	NG NUMBERS		
Rapid Screening of Meta Fluorescence Spectromet		ble High-Resolut	ion X-Ray				
6. AUTHORS				-			
Alan D. Hewitt							
7. PERFORMING ORGANIZATION NA	ME(S) AND ADDRE	ESS(ES)		4	ORMING ORGANIZATION		
U.S. Army Cold Regions 72 Lyme Road Hanover, New Hampshi		Engineering Lab	oratory		вт numвея ial Report 95-14		
9. SPONSORING/MONITORING AGEN		ADDRESS(ES)		1	NSORING/MONITORING NCY REPORT NUMBER		
U.S. Army Environmental Center Aberdeen Proving Ground, Maryland 21010-5401					SFIM-AEC-ET-CR-95039		
11. SUPPLEMENTARY NOTES							
12a. DISTRIBUTION/AVAILABILITY ST	ATEMENT			12b. DIST	RIBUTION CODE		
Approved for public rele	ase; distributi	on is unlimited.					
Available from NTIS, Spr	ingfield, Virg	inia 22161					
Analysis of copper, zinc, a cadmium, tin, and barium resolution X-ray fluoresc parameter analysis and a incoherent backscatter pe a few reference materials capability of these two rap field samples, and labora within 50% of the expected	n was perfornence spectrona second methak for matrix and are relationd XRF analytory spiked so	ned on soils and neters (XRF). Quand that relies or normalization. The ively insensitive the iss methods by depils. With the exceptions.	other particle matrice intitative determination analyte response facese two methods of into sample matrix competermining metal conception of nickel, cobal	es using tons were ctors and astrumer position tentration and characteristics.	two field-portable high- e based on fundamental d uses the Compton K_{α} atal analysis require only This study assessed the as in reference materials, aromium, concentrations		
14. SUBJECT TERMS					15. NUMBER OF PAGES 21		
Analysis Field an	alysis	Soil analysis	X-ray spectro	metry	16. PRICE CODE		
17. SECURITY CLASSIFICATION OF REPORT	18. SECURITY C	CLASSIFICATION GE	19. SECURITY CLASSIFICA OF ABSTRACT	TION	20. LIMITATION OF ABSTRACT		

UNCLASSIFIED

UNCLASSIFIED

UNCLASSIFIED

UL